


# Disposable EEG Cup+ Electrode Silver/Silver Chloride

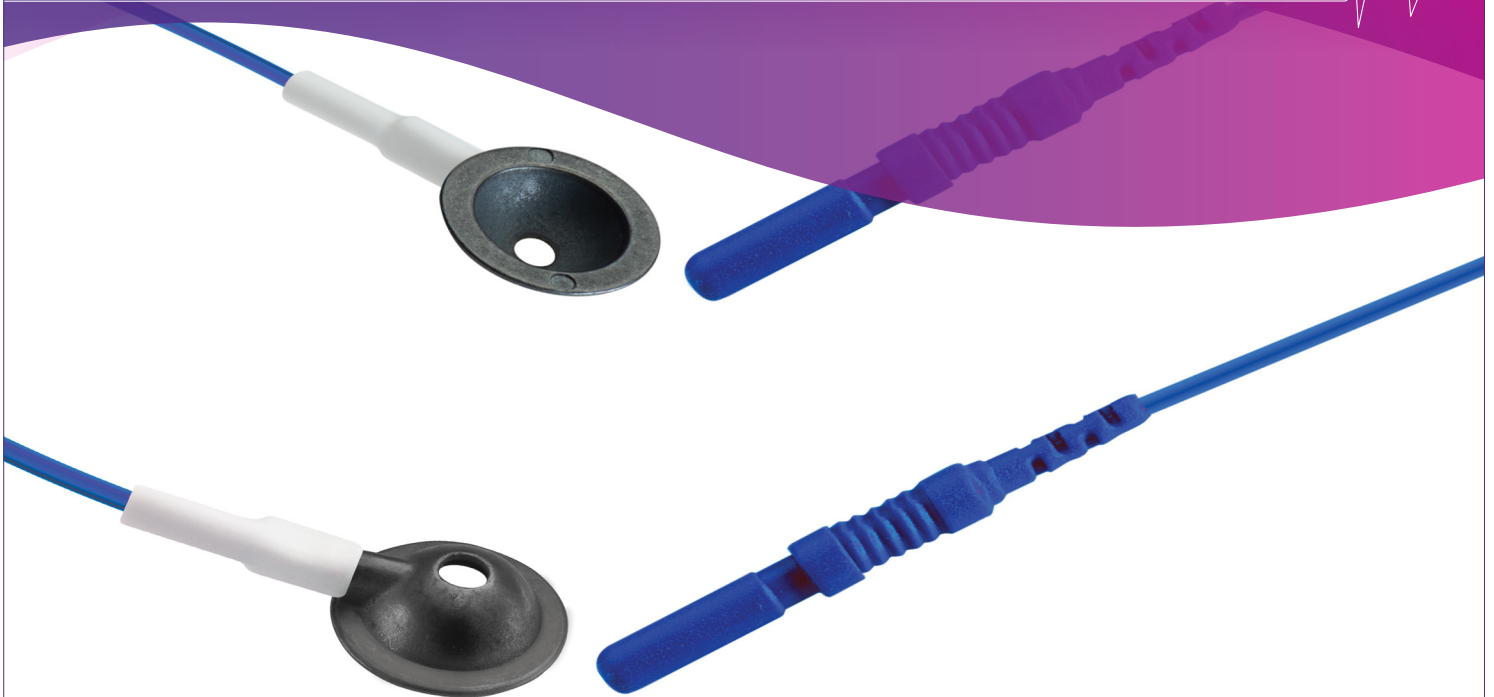

**Technomed understands the importance of clear signals, efficiency and ease of use.**

Technomed's Disposable Ag/AgCl EEG Cup+ Electrode with a larger dome holds a higher volume of conductive material.

Cup optimized for EEG procedures to focus on the patient and the recording.

The electrode is designed with flexible wires and connector optimized for easy headbox insertion and extraction.

## **User Benefits**

- No cleaning or preparation before use of the EEG electrodes needed.
- Clear and reliable signal due to optimized layer of frequency stable Ag/AgCl.
- Optimized signal quality due to larger electrode surface and volume.
- Short cup-to-cable connection to reduce lift risk.
- Optimized connector for easy headbox insertion and extraction.
- Low profile disk with soft and flexible anti-tangle lead wires for high patient comfort.
- Single use to avoid cross-patient contamination.


# Available variants

Materials	Cup Size	Cup Height	Wire Length	Packing	Part Number
Silver/Silver Chloride (Ag/AgCl)	10 mm	3 mm	1.0 m	Sets of 25, 5 x 5 colors	TE/C53-438
Silver/Silver Chloride (Ag/AgCl)	10 mm	3 mm	1.5 m	Sets of 25, 5 x 5 colors	TE/C53-638
Silver/Silver Chloride (Ag/AgCl)	10 mm	3 mm	2.0 m	Sets of 25, 5 x 5 colors	TE/C53-838
Silver/Silver Chloride (Ag/AgCl)	10 mm	3 mm	2.5 m	Sets of 25, 5 x 5 colors	TE/C53-938
Silver/Silver Chloride (Ag/AgCl)	10 mm	3 mm	1.0 m	Sets of 10, 1 x 10 colors	TE/C53-432
Silver/Silver Chloride (Ag/AgCl)	10 mm	3 mm	1.5 m	Sets of 10, 1 x 10 colors	TE/C53-632
Silver/Silver Chloride (Ag/AgCl)	10 mm	3 mm	2.0 m	Sets of 10, 1 x 10 colors	TE/C53-832
Silver/Silver Chloride (Ag/AgCl)	10 mm	3 mm	1.5 m	Sets of 15, 12 colors	TE/C53-6315
Silver/Silver Chloride (Ag/AgCl)	10 mm	3 mm	2.0 m	Sets of 15, 12 colors	TE/C53-8315
Silver/Silver Chloride (Ag/AgCl)	10 mm	3 mm	2.5 m	Sets of 15, 12 colors	TE/C53-9315

**5 colors:** red, green, white, yellow, blue

**10 colors:** red, green, white, yellow, blue, grey, purple, orange, brown, black

**12 colors:** red, 2x green, 2x white, yellow, blue, grey, purple, orange, brown, 2x black, light blue, pink

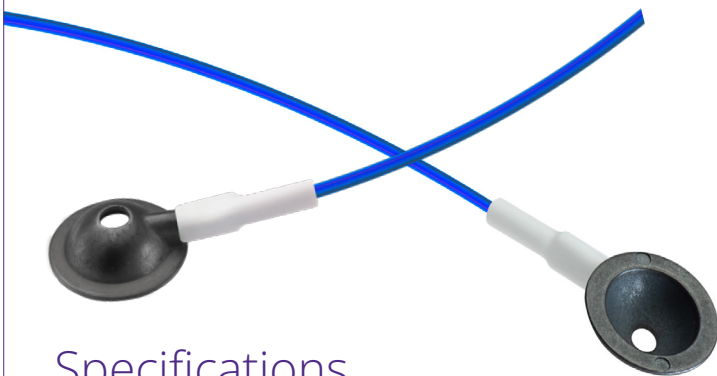

## Specifications

Cup	Ag/AgCl layer
Not made with natural rubber latex	√
REACH & RoHS (EU) compliant	√
Lead wire core	Tin Plated Copper
Lead wire insulation	PVC
Connector	DIN 42802-1

Your Technomed Business Partner:

USA  
Rx ONLY



Legal manufacturer:

Technomed Europe

Wiebachstraat 25a

6466 NG Kerkrade

The Netherlands

T +31 43 408 68 68 | infotechnomed@welcony.com

Technomed USA

1626 Terrace Drive

Roseville | MN 55113 | United States

T +1 651 288 0890 | infotechnomed@welcony.com

www.technomedmedical.com